


Manual

Please read manual before
using the product.



ELECTRO ACUPUNCTURE DEVICE

Wahonka Medical Products / CE / ISO 13485 / ISO 9001

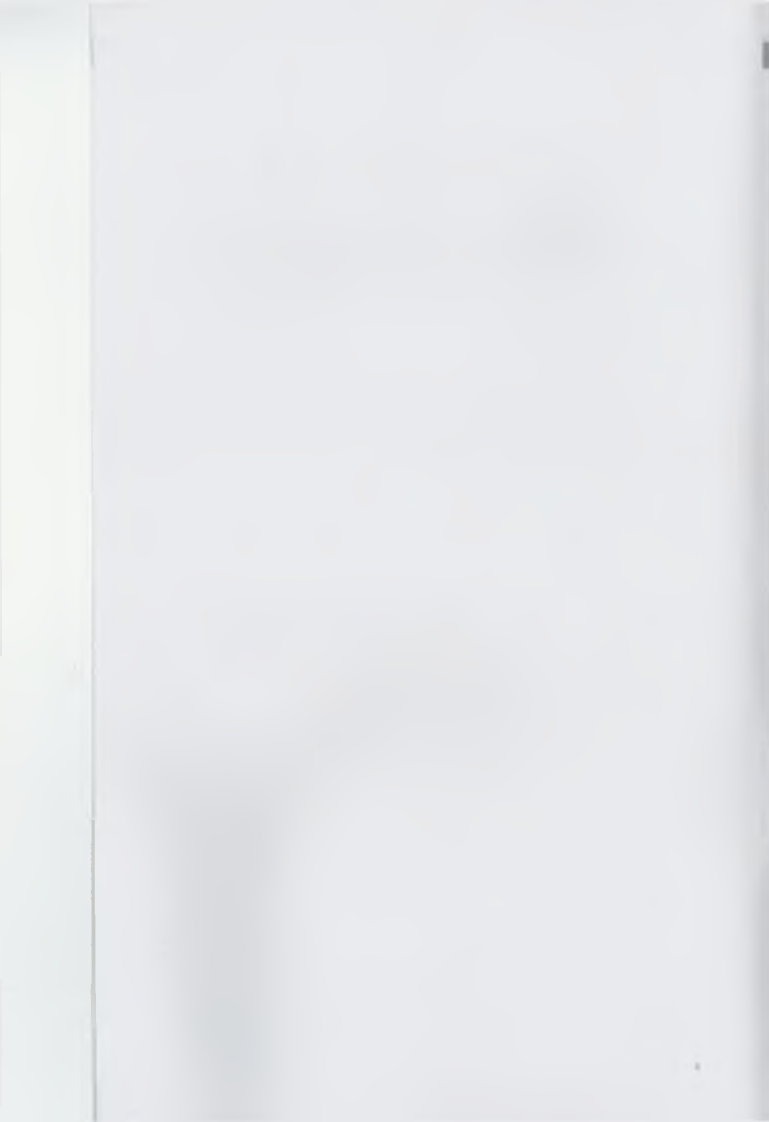

Product Introduction

Theory

Electric Acupuncture therapy device is a modern and convenient Health Device. It allows you to instantly tap into the benefits that are derived from a 5000 years old Chinese Health Science. Electric Acupuncture therapy device based on the theory that there is an invisible force or energy in nature that the Chinese call qi (pronounced as chee). In a manner similar to blood circulation, qi flows to each cell via a complex network of pathways usually called `Channels` or Meridians`. If qi gets stuck it causes an imbalance in the body. Acu points are those spots where qi has the greatest tendency to get stuck; hence stimulating them with finger pressure (Acupressure) helps to free the qi and restore the normal circulation. Once circulation of qi is restored, the cells eventually return to their normal & healthy state.

Electric Acupuncture therapy device is with digital display for location of acupuncture or trigger points. It incorporates an effective push button stimulation feature for immediate electric acupuncture treat-ment of the point or area immediately following its detection. The digital display helps the user to locate the most accurate acupuncture or trigger point easily, as the number increases with the lowering of the electric resistance of the skin (when a maximum number is reached, an acupuncture point is detected and electroacupressure may begin). Self-treatment for electric acupuncture has never been so easy!

Electric Acupuncture therapy device has additional features, such as adjustable frequency from 2 - 23Hz,

It is one of the most innovative hand-held stimulators available in the market. When the metal probe touches a low resistance point the loudspeaker gives a sound signal and the display screen show the number, using the fingers to skillfully press the key points, which stimulate the body's natural self-curative abilities. When these trigger points are pressed they release muscular tension, promote circulation of blood and the body's life force energy to aid healing. Acupressure uses gentle to firm pressure and integrates bodywork therapies, therapeutic touch healing, imagery, energy psychology and therapeutic massage therapy techniques.

Specification

Output Intensity	0 to 10 mA
Output Intensity Control	ON / OFF, 0 to 10
Power Supply	DC 3.7V
Micro-current Output	0 to 2 microamps
Sensitivity Adjust Control	0 to 10
Frequency Control	2-23Hz(Adjustable)
Large LCD Digital Display	Frequency/ Sensitivity Indicator
Indicator Lamp	Pulse Rate Indicator
Pulse Width	220 microseconds
Pulse Modes	Continuous
Pulse Shape	Asymmetric biphasic square wave
Point Detection	Meter, Sound and Light Indicators
Output Channels	One
Box Size	269 x 170 x 53 mm

Replaceable probe
of finding acupoints

power and
intensity control

display screen

treatment button

sensitive
control

light



Ashi pads
interface

metal
conduction 1

metal
conduction 2

frequency
control



Conductive stick
interface

Charging Interface



Cable

Beauty probe

Mini probe

Replaceable probe



Conductive bar

Product features

1. Automatically find points, intelligent and easy.
2. Good for blood circulation.
3. portable and easy to use.
4. Adjusting the sensitivity of points detection according to your skin type and skin condition.
5. The user-friendly design, easy to operate.
6. Audio Tips.
7. We use lithium batteries that can be recharged.
8. Treatment probe can be replaced.

Scope of application


Suitable for human wrist, elbow, ankle, knee, neck, waist, face and other parts have points.

Warning

1. Not using with electronic devices, artificial heart, lung and other medical electronic equipment at the same time.
2. Pregnant women and women with acute menstrual period, severe heart disease, severe abnormal blood pressure, patients with severe hypertension, skin abnormalities, fever and other prohibited
3. Do not give children or people without the ability to use the expression of consciousness
4. cancer and other severe, prohibit the use of emergency patients.

Unique Features

1. Adjusting the "switch / intensity" button to 4, and then adjusting the "sensitivity" button to about 2, the detection sounds 2 times per second, the frequency can be adjusted between 2-23Hz, you can adjust the intensity according to your own feelings. Because everyone's skin is different and the humidity of weather is different, "sensitivity" adjusting vary.

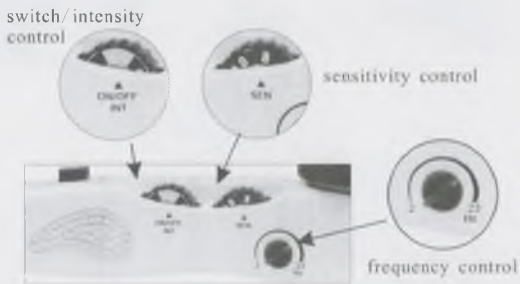

2. Before using, the treatment probe is needed to be mounted on a host, and then insert the conductive wire rod to the host, the other hand holding the stick to form a loop.




3. Hold the acupuncture pen with a correct posture, like holding a pen, touching the metal conductive sheet 1 with your thumb, and touching the metal conductive sheet 2 with your Hukou acupoint (please see the picture below), and then put the probe to touch your skin uprightly.



4. When probing the meridian points and when you hear a beep and display value fast connections increases, the meridian points are detected, then press the button to start the therapy treatment. Note: When detecting the meridian points, can not press the treatment button.



5. When doing the treatment, you can adjust the "switch / intensity" button to the appropriate treatment intensity, to feel comfortable and appropriate (a total of 10 former level is adjustable from the digital 0-10).

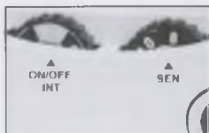



6. After using, turn the "switch / intensity" button to "off".



7. When replacing the probe to beauty probe, you need to put off the former probe anticlockwise and put on the beauty probe clockwise.



8. When charging, the LCD will appear "CH6" words and flashing green light, when fully charged, the green light is off.

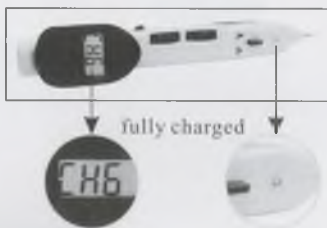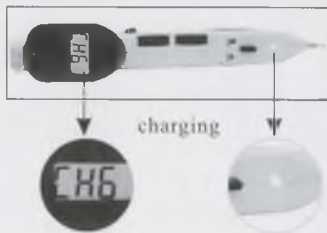

Treatment time

One or twice a day, For treatment in one acupoint or one area of skin is around 3-5minutes.

Cleaning and maintenance

About cleaning:

To wipe the surface with a wet soft towel.

Probe maintenance:

- 1, When the probe's surface gets dirt, you can wipe the probe with a towel.
- 2, In order to avoid cross-infection, when other people finish using, it is recommended to use alcohol to wipe the probe.